

Corrosive Cabinet

WA810120B



It has acquired FM & CE Certification and is in strict comply with FM Approval Standard for Storage Cabinets (Flammable and Combustible Liquids) Class Number 6050, OSHA 29 CER 1910.106 of American Occupational Safety and Health Administration, NFPA CODE30 of American National Fire Protection Association and Hazardous Chemicals Control Ordinance of State Bureau of Safe Production Supervision and Administration

SYSBEL Corrosive Cabinets for storing corrosive liquids, Apply to the LABS of colleges and universities, research institutes, etc.

- 1.The standard PVC tray to prevent leakage of weak corrosive liquid corrosion cabinet put oneself in another's position.
- 2.Specially designed small size safety cabinet easy to any place on the stage or the audience (laboratory).
- 3.Double-wall construction with 38mm insulating air space for fire resistance.
- 4.Over 1.2mm thick, fully welded,construction holds squareness for longer life, offering greater protection in a fire.
- 5.5.3cm leak tight sump at the bottom of cabinet's maximally catches incidental drips.
- 6.Door can be fully opened to 180°, easy to operate, three-point latch with a manual lock for better security.
- 7.Standardized warning label is highly visible and anti-corrosive.
- 8.Unique spill-catcher shelves catch incidental drips and adjust on 6cm centers.
- 9.Durable&chemical resistant,lead-free powder coat inside and outside cabinet,reduce effects of corrosion and humidity.
- 10.2inches vents with integral flame arresters in both sides of every cabinet.
- 11.In terms of OSHA, on outside side panel,there is built-in grounding static connector for easy grounding.



China's first certified by FM




Corrosive Cabinet WA810120B

Volume (Gal/L)	Certification	Door type	Ext dimension (H×W×D/cm)	Packing dimension (H×W×D/cm)	Adjustable Shelves(Pcs)	Shelf Loading Weight(Kgs.)	N.W. (Kgs.)	G.W. (Kgs.)	Model
12/45	FM, CE	Single Door/ Manual	89x59x46	103x66x53	1	50	48	60	WA810120B


WA810120B



Product details

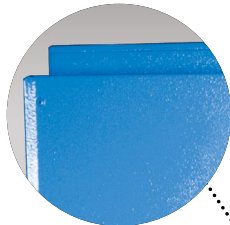

18Ga(1-mm)quality galvanized steel
Unlead coating with polish surface,
which is corrosion and moisture proof

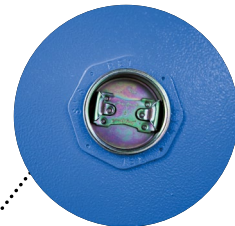

Fireproofing vents to effectively reduce
the concentration of volatile substances,
and prevent fire hazards. Special spill-
proof shelves to avoid chemical leakage
by accident.




Adjustable shelves with anti-
leakage design




The cabinet has three-point linkage lock,
and the door can open 180°
Two person management in line with
Hazardous Chemicals Control Ordinance
(December 1st, 2011 Version)



Shelf bracket welded on both
sides of cabinet

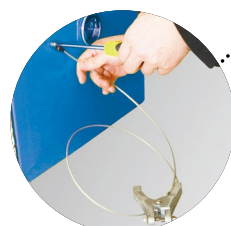

FM Approved Anti-Static
Device: Conductor grip combine
the equipment with the cabinet to
conduct the static to earth, thus to
prevent fire accident

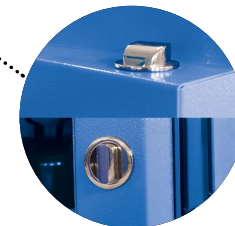

The three-point bullets: Anti-static design
structure

